

# CENTRIX Encore® D/C MiniMix®

Unit Dose Dual-Cure Core Buildup Material

## Makes Core Buildup Easier

The MiniMix unit-dose technology eliminates waste and cross-contamination concerns.



### Single-dose. Great price. Problems solved.

Our new single-dose Encore D/C MiniMix kit contains 36 cartridges. Each cartridge has two mixing nozzles and our exclusive needletip. That's 72 core buildups! Our product is dual-cure, for fast light-activation. Plus, the smaller size and ease of use will save you valuable chair time and reduce patient stress.

#### Features:

- Light-activate for 40 seconds for complete cure
- Needletip allows delivery of material to bottom of a post hole
- MiniMix technology virtually eliminates waste
- Non-slumping formula so there is no need to use a matrix
- Cuts like dentin to reduce ditching
- Radiopaque



**Encore® D/C MiniMix®  
is the only dual-cure  
sold by dose.**



**centrix®**

Making Dentistry Easier.™

CE  
0470


**6 -steps for building a void-free, bonded, composite post and core:**



1. Enlarge root canal to accept Encore D/C needletip. Bevel the edges to eliminate sharp angles.




2. Etch the post hole, rinse thoroughly and dry.

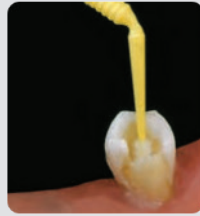

3. Apply bonding agent according to manufacturer's instructions.



4. Insert needletip to the bottom of the post hole and inject material with slow steady pressure.

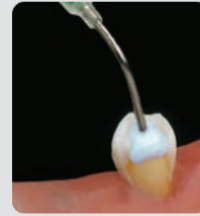

5. Continue filling by injecting and slowly withdrawing the needletip.



6. If desired, light-activate Encore D/C to accelerate the set. Allow material to fully cure and then cut and shape the core. The result is a solid, securely bonded post and core.



**Encore D/C MiniMix Specifications**

Comprehensive Strength	222 Mpa
Flexural Strength	110 Mpa
Diametral Tensile Strength	48 Mpa
Light-activated Depth of Cure	2mm 40 Seconds

**ORDERING INFORMATION**


**Encore D/C MiniMix**

Contains 36 single-dose cartridges (0.5mL net ea.) and 72 mixing nozzles.

**REF 310254** Natural Kit

**REF 310255** Contrast Blue Kit

**Shade Reference**

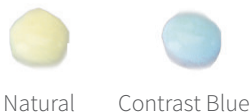

**Suggested Accessories**

**REF 130651** Snub Nose Snap-Fit Gun

**REF 910630** MiniMix Nozzles, Clear Hub, 10-pack

**Ordering is easy. Call 800.235.5862 or contact your preferred dealer. Visit [centrixdental.com](http://centrixdental.com) for more information.**